

WHY CHOOSE EUROGUARD® TRAY WRAPS?

- The SMMS material is very user-friendly, is treated with antistatic and folds very smoothly.
- The tray wrap combination 40 NAR with 50 AR is both a cost effective as practical solution, the majority of the trays can be packed with this interleaved combination.
- Excellent quality of the SMMS material, in line with the rest of our Euroguard® program.
- The blue material is alcohol repellent and is therefore also usable as sterile field.
- Possible damage of the outer wrap is immediately visible because of the green/blue combination of materials.
- As mentioned above the **EUROGUARD® TRAY WRAPS** meet all requirements of the EN-ISO 11607. A technical datasheet is available on request.
- The tray wraps can be used for steam and ETO sterilisation.

Euroguard® Tray Wrap SMMS 40 grams green

		Measuring Metric	Value	Norm EN868 Minimum
Tensile strength MD	EN-ISO 1924-2	N/5cm	70-75	55
Tensile strength CD	EN-ISO 1924-2	N/5cm	45-50	30
Tear resistance MD	EN-21974	mN	>750	750
Tear resistance CD	EN-21974	mN	>1000	1000

Euroguard® Tray Wrap SMMS 50 grams blue

		Measuring Metric	Value	Norm EN868 Minimum
Tensile strength MD	EN-ISO 1924-2	N/5cm	85-90	70
Tensile strength CD	EN-ISO 1924-2	N/5cm	45-50	30
Tear resistance MD	EN-21974	mN	>750	750
Tear resistance CD	EN-21974	mN	>1000	1000




MEDICA
EUROPE

Galliërsweg 20, 5349 AT Oss +
P.O. Box 746, 5340 AS Oss +
The Netherlands +
T +31 412 671300 +
F +31 412 671310 +
www.medica-europe.nl +

Versie 3 | 10-2019

EUROGUARD®

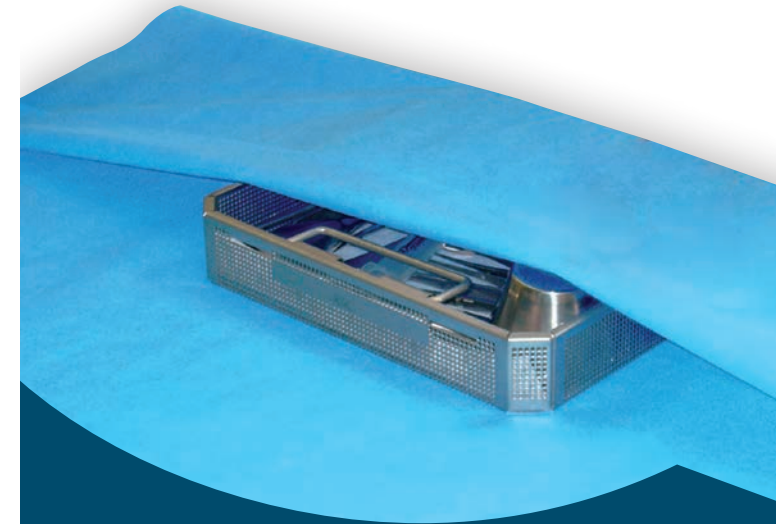

EUROGUARD® TRAY WRAPS



TRAY WRAPS

Medica Europe develops, manufactures and markets a wide range of medical disposables for the Operating Room, cardiac catheterisation laboratories, critical care units and interventional radiology centers. Since many years we manufacture drapes and gowns under our own brand Euroguard®. We have become specialists in the SMS and SMMS materials. With this expert knowledge we extended our range to the CSR field and we developed a sterilization packaging under the name **EUROGUARD® TRAY WRAPS**.

Medica Europe is aware of the responsibilities of the CSSD and delivers only sterilization packaging which meets the European standards. The **EUROGUARD® TRAY WRAPS** meet the EN-ISO 11607.

EUROGUARD® TRAY WRAPS are made of SMMS material and consist of four layers of polypropylene (Spunbound, Meltblown, Meltblown, Spunbound). The Spunbound layers are melted and spun, the fibres give the strength but also the drapeability characteristics. The Meltblown layers consists of microscopic little fibres which create a web and frame the barrier for micro organism.

EUROGUARD® TRAY WRAPS is available in one interleaved combination:

- 40 grams green, non alcohol repellent, outer wrap
- 50 grams blue, alcohol repellent, inner wrap

EUROGUARD® TRAY WRAP blue is resistant to alcohol. The SMMS material has excellent moisture repellency which prevents the migration of micro organism.

The tray wraps are available in different sizes to fit the different dimensions of the trays:

75 x 75 cm

90 x 90 cm

100 x 100 cm

120 x 120 cm

